

HiRes Agarose

PUREGENE

Product code : PG-432

An agarose for molecular screening that improves resolution of small DNA fragments and PCR products.

The key to producing the HiRes (High Resolution) agaroses is harvesting the appropriate seaweed at precise time in its growth cycle. There are also certain modifications in the chemical structure of the polymer during the manufacturing process. Puregene has come up with HiRes Agarose for applications that require efficient separation of small DNA fragments and PCR products.

FEATURES

High resolution of short PCR products and DNA fragments.

Improved clarity of the gel, enhancing visibility.

Better handling than competitive products because of a stronger gel structure and higher gel strength. The chances of gels breaking or cracking when handled are greatly minimized, even with lower concentrations of agarose.

High gel strength allows use in blotting.


FUNCTIONAL TESTS

DNA resolution: bands appear sharp and finely resolved.

DNase/RNase activity: none detected.

Gel background: very low after EtBr staining.

DNA binding: very low.

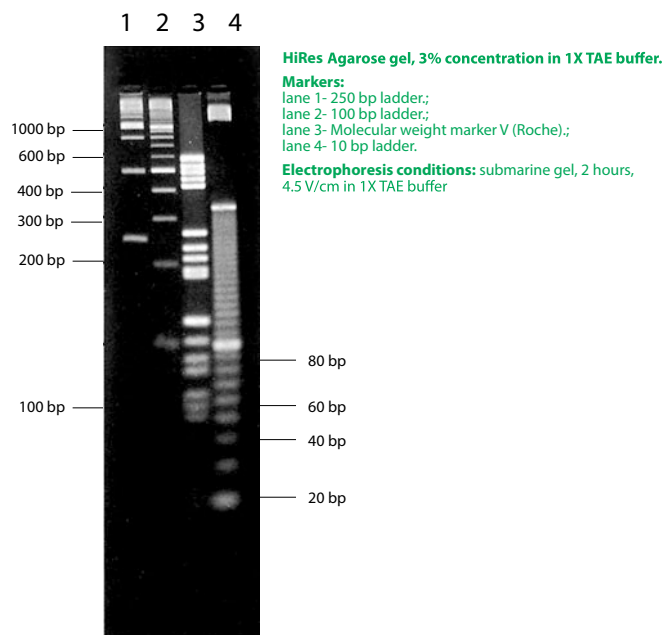

SPECIFICATIONS

*EEO((electroendosmosis)

	1.5%	3%
Moisture	≤ 10%	
Ash	≤ 0.35%	
EEO*	≤ 0.12%	
Sulfate	≤ 0.12%	
Clarity (NTU)	≤ 5	
Gel Strength (g/cm ²)	≥ 600	≥ 1500
Gelling Temperature (°C)	≤ 35.5	
Melting Temperature (°C)		

RANGES OF SEPARATION

1.8%	400 – 1200 bp
3.0%	150 – 800 bp
4.5%	15 – 400 bp

These ranges are approximate and have been calculated in 1X TAE buffer. To achieve the best resolution of HiRes agarose gels, they should be stored at 4° - 8° C for 30 minutes before use.