

INTEGRA



PIPETBOY acu Light, coloured and always ready to use

PIPETBOY acu

Light, coloured and always ready to use



Colour-coded

You can select your own favourite colour for your personal PIPETBOY acu and bring some colour into your lab! In this way you can work more safely and avoid mix-ups.

You can allocate colours to particular applications. For example: dedicated instruments for particular cell lines, virology or aspiration of waste products.



Total Versatility

PIPETBOY acu is ideal for use with all commercially available serological pipettes. The pipette mount made of silicone provides a safe and secure grip. It is durable, flexible and autoclavable.



Precise dosing

Together with the speed control, the dosing knob guarantees precise, fine dosing.



Battery indicator

An illuminated red LED shows when the battery needs to be charged.



green



red



blue



transparent



purple



Light as a feather

The PIPETBOY acu is light as a feather and ergonomically designed for ease of handling.



Variable pipetting speed

Control the pipetting speed continuously via the knurled screw.



Always on standby

No more waiting for charge-up. This can be done whilst you are using the instrument. The high performance and environmentally friendly NiMH battery technology incorporated in PIPETBOY acu makes it constantly available for use – and the battery shows no "memory effect".



Maximum safety

The easy-to-change sterile filter provides your PIPETBOY acu with reliable protection from liquid contamination and damage. And if you do happen to aspirate critical liquids with the pipette, the silicone pipette mount, the pipettor nose piece and the filter rubber are autoclavable while the sterile filter can be replaced. So easy and so safe.

PIPETBOY acu

Ordering information

Description	Part No.
PIPETBOY acu classic , incl. wall mount, mains adapter, sterile filter, 9-V NiMH-battery	155 000
PIPETBOY acu green , as Art. No. 155 000	155 015
PIPETBOY acu red , as Art. No. 155 000	155 016
PIPETBOY acu blue , as Art. No. 155 000	155 017
PIPETBOY acu transparent , as Art. No. 155 000	155 018
PIPETBOY acu purple , as Art. No. 155 000	155 019

Ordering information for accessories


No.	Description	Part No.
1	Battery compartment lid for PIPETBOY acu classic	155203
	Battery compartment lid for PIPETBOY acu green	155160
	Battery compartment lid for PIPETBOY acu red	155161
	Battery compartment lid for PIPETBOY acu blue	155162
	Battery compartment lid for PIPETBOY acu transparent	155163
	Battery compartment lid for PIPETBOY acu purple	155164
2	O-ring	153 235
	Pipette mount, grey, incl. O-ring for PIPETBOY acu green, blue, red and purple	155 060
3	Pipette mount, green, incl. O-ring for PIPETBOY acu classic and transparent	155 230
	Filter rubber	153 225
4	Filter sterile (pore size 0,45 µm)	153 015
	Filter unsterile (pore size 0,45 µm)	153 016
	Filter sterile (pore size 0,2 µm)	156 608
	Filter unsterile (pore size 0,2 µm)	156 607
	Pipette mount rubber	151 020
7	Sterile module, complete for PIPETBOY acu classic and transparent	155 025
	Sterile module, complete for PIPETBOY acu green, blue, red and purple	155 070
	9 V NiMH-battery	153 005
8	Separate rapid charger for two 9 V batteries (including one battery)	153 010
	Mains Adapter	EU version 153 210
	(100 – 240 VAC, 50/60 Hz):	US version 153 211
		UK version 153 214
	Wall Mount	155 521
9	Stand for PIPETBOY (acu or comfort) and VACUBOY	155 065

Technical data

Dimensions (LxHxW):	125 x 130 x 35 mm
Weight:	220 g
Mains adapter output:	17 VDC/180 mA
Batteries:	1 x 9 V NiMH 180 mAh
Operating condition:	5 – 40° C, max. 80% RH



The PIPETBOY acu is constructed in such a way that parts can be replaced on the spot. This is a rapid procedure, making PIPETBOY acu ready for work immediately. The spare parts illustrated in the graphic can be ordered individually.

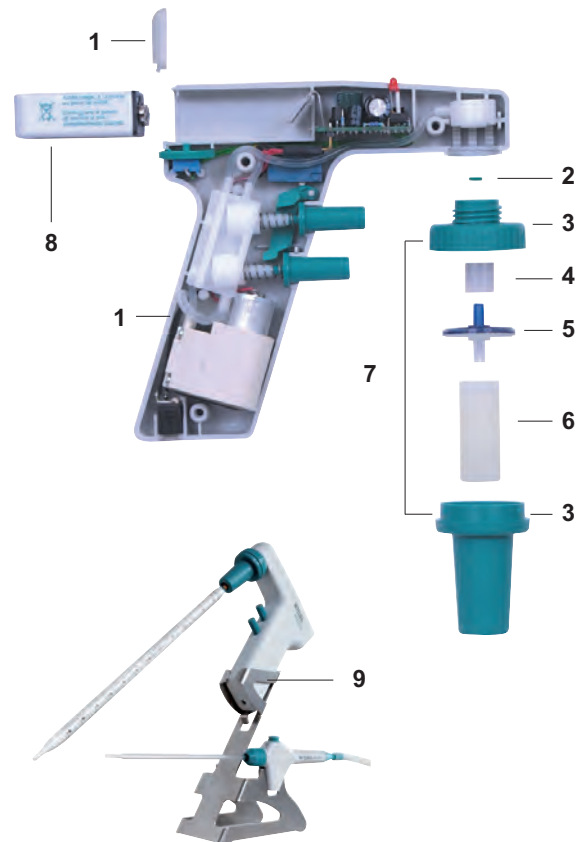

CE
Swiss made

INTEGRA Biosciences AG

CH-7205 Zizers, Switzerland
Phone: +41 81 286 95 30
Fax: +41 81 286 95 33
info@integra-biosciences.com
www.integra-biosciences.com



INTEGRA