

PUREGENE™
GENETIX BRAND

4 COLOR PRESTAINED PROTEIN LADDER



Cat # PG-PMT0782

Size: 2X250µl

Store at -20°C

PUREGENE™
GENETIX BRAND

FEATURES:

Wide Range: With 10 recombinant proteins covering ~10-180kDa

Ready to Use: Supplied in a loading buffer for direct loading on gels.

Easy to Identify: One pink color reference bands at ~137kDa, an orange color at ~71kDa and one green reference band at ~51kDa.

(Lot to lot variation of the apparent molecular weight of prestained protein ladder is ~5%)


Sharp Bands

APPLICATIONS:

- Monitoring of protein separation during SDS-PAGE.
- Verifying Western transfer efficiency.
- Approximate sizing of proteins on SDS-polyacrylamide gels and Western blots.
- Locating a protein of interest for excision from an unstained preparative gel.

QUALITY CONTROL:

Tested in SDS-Polyacrylamide Gel electrophoresis and Western blotting.

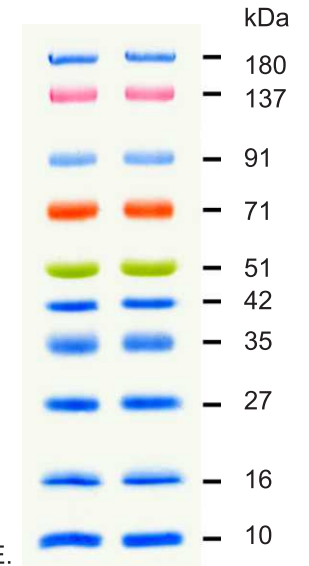

7µl of the Ladder is analyzed on 4-20% SDS PAGE gel in Tris-Glycine Buffer

Storage

Stable at +4°C for 3 months
Store at -20°C for 12 months

 **genetix**
Biotech Asia Pvt. Ltd.

MFG

71/1, First Floor, Shivaji Marg, New Delhi - 110015 INDIA
Phones : +91-11-45027000, Fax : +91-11-2541-9631
E-mail : info@genetixbiotech.com Website: www.genetixbiotech.com